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Εναρμονισμένα Πρότυπα	Σύντομη περιγραφή
Common specifications	Non available for the product
DIN EN ISO 13485	Medical devices - Quality management systems - Requirements for regulatory purposes
DIN EN 62366-1	Medical devices - Part 1: Application of usability engineering to medical devices
DIN EN ISO 14971	Medical devices – Application of risk management to medical devices
DIN EN ISO 14155	Clinical Investigation of Medical Devices for Human Subjects - Good Clinical Practice
DIN EN ISO 10993-1	Biological evaluation of medical devices - Part 1: Evaluation and testing within a risk management system
DIN EN ISO 10993-3	Biological evaluation of medical devices - Part 3: Tests for genotoxicity, carcinogenicity and reproductive toxicity
DIN EN ISO 10993-5	Biological evaluation of medical devices - Part 5: Tests for in vitro cytotoxicity
DIN EN ISO 10993-6	Biological evaluation of medical devices – Part 6: Tests for local effects after implantation
DIN EN ISO 10993-10	Biological evaluation of medical devices - Part 10: Tests for irritation and skin sensitization
DIN EN ISO 10993-11	Biological evaluation of medical devices - Part 11: Tests for systemic toxicity
DIN EN ISO 10993-12	Biological evaluation of medical devices - Part 12: Sample preparation and reference materials
DIN EN ISO 11737-1	Sterilization of medical devices - Requirements for the estimation of population of microorganisms on a product
DIN EN ISO 11737-2	Sterilization of medical devices - Microbiological methods - Part 2: Tests of sterility performed in the definition, validation and maintenance of a sterilization process
DIN EN ISO 11137-1	Sterilization of health care products - Radiation - Part 1: Requirements for development, validation and routine control of a sterilization process for medical devices
DIN EN ISO 11137-2	Sterilization of Health Care Products - Radiation - Part 2: Establishing The Sterilization Dose
DIN EN 556-1	Sterilization of medical devices - Requirements for medical devices to be designated "STERILE" - Part 1: Requirements for terminally sterilized medical devices
DIN EN ISO 11607-1	Packaging for terminally sterilized medical devices - Part 1: Requirements for materials, sterile barrier systems and packaging systems
DIN EN ISO 11607-2	Packaging for terminally sterilized medical devices - Part 2: Validation requirements for forming, sealing and assembly processes
DIN EN ISO 20417	Medical devices – Information to be supplied by the manufacturer
DIN EN ISO 15223-1	Medical Devices - Symbols To Be Used With Medical Device Labels, Labelling And Information To Be Supplied - Part 1: General Requirements
DIN EN 868-2	Packaging for terminally sterilized medical devices - Part 2: Sterilization wrap - Requirements and test methods
DIN EN 868-5	Packaging for terminally sterilized medical devices - Part 5: Sealable pouches and reels of porous materials and plastic film construction - Requirements and test methods
DIN EN ISO 14698-1	Cleanrooms and associated controlled environments -- Biocontamination control -- Part 1: General principles and methods

DIN EN ISO 14698-2	Cleanrooms and associated controlled environments -- Biocontamination control -- Part 2: Evaluation and interpretation of biocontamination data
ISTA 2a	Partial Simulation Performance Tests - Packaged Products 150 lb (68 kg) or Less
USP <151>	Pyrogen Study
ASTM F1886/F1886M	Standard Test Method for Determining Integrity of Seals for Flexible Packaging by Visual Inspection
ASTM F88/F88M	Standard Test Method for Seal Strength of Flexible Barrier Materials
ASTM F3039	Standard Test Method for Detecting Leaks in Nonporous Packaging or Flexible Barrier Materials by Dye Penetration

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Albertsson, A., Eklund, M. (1995). Influence of Molecular Structure on the degradation mechanism of Degradable Polymers: In Vitro Degradation of Poly(Trimethylene Carbonate, Poly(Trimethylene Carbonate-co-Caprolactone), and Poly(Adipic Anhydride). *J Applied Polymer Sci* 57, 87-103.

Ali, S. A., Zhong, S. P., Doherty, P. J. & Williams, D. F. (1993). Mechanisms of polymer degradation in implantable devices. I. Poly(caprolactone). *Biomaterials*, 14, 648-56.

Apelqvist, J., Bakker, K., Van Houtum, W. H., Nabuurs-Franssen, M. H. & Schaper, N. C. (2000). International consensus and practical guidelines on the management and the prevention of the diabetic foot. International Working Group on the Diabetic Foot. *Diabetes Metab Res Rev*, 16 Suppl 1, S84-92.

Baartmans, M. G., Dokter, J., Den Hollander, J. C., Kroon, A. A. & Oranje, A. P. (2011). Use of skin substitute dressings in the treatment of staphylococcal scalded skin syndrome in neonates and young infants. *Neonatology*, 100, 9-13.


Barbachowska A, Korzeniowski T, Surowiecka A, Tomaka P, Bugaj-Tobiasz M, łączyk M, Górecka Z, Chrapusta A, Strużyna J. The Effectiveness of an Alloplastic Epidermal Substitute in the Treatment of Burn Wounds in Children: A Comparative Clinical Study of Skin Substitutes and Silver and Paraffin Gauze Dressings. *J Clin Med*. 2024 Nov 28;13(23):7238.

Blome-Eberwein, S.A., Amani, H., Lozano, D.D., Gogal, C., Boorse, D., Pagella, P. Burns. A bio-degradable synthetic membrane to treat superficial and deep second degree burn wounds in adults and children – 4 year experience. 2020 Aug 29;S0305-4179(20)30507-6

Bostman, O. M. (1991). Absorbable implants for the fixation of fractures. *J Bone Joint Surg Am*, 73, 148-53.

Brady, J. M., Cutright, D. E., Miller, R. A. & Barristone, G. C. (1973). Resorption rate, route, route of elimination, and ultrastructure of the implant site of polylactic acid in the abdominal wall of the rat. *J Biomed Mater Res*, 7, 155-66.

Bryant, R., Nix, D. (2006). *Acute and Chronic Wounds* 3rd Edition.

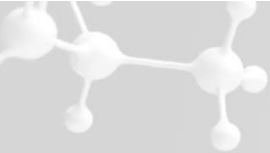Burd, A. & Yuen, C. (2005). A global study of hospitalized paediatric burn patients. *Burns*, 31, 432-8.

Cha, Y. & Pitt, C. G. (1990). The biodegradability of polyester blends. *Biomaterials*, 11, 108-12.

Cheema L, Manzoor S, Khalid U, Shamim R, Hashaam, Tayyab Z, Bashir M. Suprathel Dressing at Split Thickness Skin Graft Donor Site for Pain Control and Wound Healing. *Pakistan journal of medical and health sciences*, 2022, 16(10), 116-118 Dastagir N, Kijas D, Obed D, Tamulevicius M, Vogt PM, Dastagir K. Suprathel® and water-filtered infrared-A radiation (wIRA) as a new treatment strategy for toxic epidermal necrolysis (TEN): A prospective study. *Burns*. 2024 Dec;50(9):107283

Delgado-Miguel C, García Morán A, Fuentes Gómez L, Díaz M, Miguel-Ferrero M, López-Gutiérrez JC. Comparison of the effectiveness of three different skin substitutes for the treatment of pediatric burns. *Eur J Pediatr*. 2024 Dec 13;184(1):80.

Demidova O, Manushin S. Alloplastic skin substitute SUPRATHEL® dressings in treatment of donor sites in children with burns. EBA 2017, Barcelona, Spain (Available from: <https://www.morressier.com/article/alloplastic-skin-substitute-dressings-treatment-donor-sites-children-burns/594bbebfd462b8028d893e61>)

Demircan, M., Gürünlüoğlu, K., Gözükar, Bag, H.G., Koçbıyık, A., Gül, M., Üremiş, N., Gül, S., Gürünlüoğlu, S., Türköz, Y., Taşçı, A. Impaction of the polylactic membrane or hydrofiber with silver dressings on the interleukin-6, tumor necrosis factor- α , transforming growth factor- β 3 levels in the blood and tissues of pediatric patients with burns. *Ulus Travma Acil Cerrahi Derg*. 2021 Jan;27(1):122-131

Dhivya, S., Padma, V. V. & Santhini, E. (2015). Wound dressings - a review. *Biomedicine (Taipei)*, 5, 22.

Everett, M., Massand, S., Davis, W., Burkey, B. & Glat, P. M. (2015). Use of a copolymer dressing on superficial and partial-thickness burns in a paediatric population. *J Wound Care*, 24, S4-8.

Fernandes S, Teixeira I, Carmo L, Campos M, Garcia M. The Use of a Polylactic Membrane in Pediatric Burns Proves to be Successful Even After Late Application. *J Burn Care Res*. . 2023 Sep 7;44(5):1176-1181

Fischer, S., Kremer, T., Horter, J., Schaefer, A., Ziegler, B., Kneser, U., Hirche, C. Suprathel® for severe burns in the elderly: Case report and review of the literature. *Burns*. 2016 Aug;42(5):e86-92


Galati V, Vonthein R, Stang F, Mailaender P, Kisch T. Split thickness skin graft versus application of the temporary skin substitute suprathel in the treatment of deep dermal hand burns: a retrospective cohort study of scar elasticity and perfusion. *Int J Burns Trauma*. 2021 Aug 15;11(4):312-320

Gurunluoglu, K., Demircan, M., Tasci, A., Uremis, M. M., Turkoz, Y., Bag, H. G., Akinci, A. & Bayrakci, E. (2019). The Effects of Two Different Burn Dressings on Serum Oxidative Stress Indicators in Children with Partial Burn. *J Burn Care Res*, 40, 444-450.

Gürünlüoğlu K, Demircan M, Koç A, Koçbıyık A, Taşçı A, Durmuş K, Gürünlüoğlu S, Gözükar Bağ H. The Effects of Different Burn Dressings on Length of Telomere and Expression of Telomerase in Children With Thermal Burns. *J Burn Care Res*. 2019 Apr 26;40(3):302-311.

Hakkarainen, T., Koivuniemi, R., Kosonen, M., Escobedo-Lucea, C., Sanz-Garcia, A., Vuola, J., Valtonen, J., Tammela, P., Mäkitie, A., Luukko, K., Yliperttula, M., Kavola, H. J Control Release. Nanofibrillar cellulose wound dressing in skin graft donor site treatment. 2016 Dec 28;244(Pt B):292-301.

Harenberg, P. S., Hrabowski, M., Ryssel, H., Gazyakan, E., Germann, G., Engel, H. & Reichenberger, M. A. (2010). CASE REPORT Febrile Ulceronecrotic Mucha-Habermann Disease. *Eplasty*, 10.

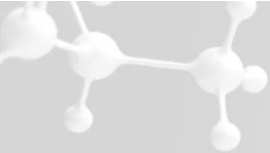Summary of Safety and Clinical Performance (SSCP)

based on Clinical
Evaluation Report V11

issued: 2025.02.28

Heitzmann W, Mossing M, Fuchs PC, Akkan J, Seyhan H, Grieb G, Opländer C, Schiefer JL. Comparative Clinical Study of Suprathel® and Jelonet® Wound Dressings in Burn Wound Healing after Enzymatic Debridement.

Hettiaratchy, S. & Papini, R. (2004). Initial management of a major burn: II--assessment and resuscitation. *BMJ*, 329, 101-3.

Highton, L., Wallace, C., Shah, M. Use of Suprathel® for partial thickness burns in children. *Burns*. 2013 Feb;39(1):136-41

Hollinger, J. O. & Battistone, G. C. (1986). Biodegradable bone repair materials. Synthetic polymers and ceramics. *Clin Orthop Relat Res*, 290-305.

Hundeshagen, G., Collins, V. N., Wurzer, P., Sherman, W., Voigt, C. D., Cambiaso-Daniel, J., Nunez Lopez, O., Sheaffer, J., Herndon, D. N., Finnerty, C. C. & Branski, L. K. (2018). A Prospective, Randomized, Controlled Trial Comparing the Outpatient Treatment of Pediatric and Adult Partial-Thickness Burns with Suprathel or Mepilex Ag. *J Burn Care Res*, 39, 261-267.

Hunt, T., Aslam, R., Beckert, S., Wagner, S., Ghani, R., Hussain, M., Roy, S., And Sen, C. (2007). Aerobically Derived Lactate Stimulates Revascularization And Tissue Repair Via Redox Mechanisms. *Antioxidants & Redox Signaling Volume 9, Number 8*

International consensus. Wound infection in clinical practice. *Int Wound J*. 2008 Jun;5 Suppl 3(Suppl 3):iii-11. doi: 10.1111/j.1742-481X.2008.00488.x. PMID: 18489408; PMCID: PMC7951552.

Kaartinen, I. S. & Kuokkanen, H. O. (2011). Suprathel((R)) causes less bleeding and scarring than Mepilex((R)) Transfer in the treatment of donor sites of split-thickness skin grafts. *J Plast Surg Hand Surg*, 45, 200-3.

Karlsson M, Steinvall I, Elmasry M. Suprathel® or Mepilex® Ag for treatment of partial thickness burns in children: A case control study. *Burns*. 2023 Mar 11;S0305-4179(23)00043-8.

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Kamolz, L., Lumenta, D., Kitzinger, H., Frey, M. (2008). Tissue engineering for cutaneous wounds: an overview of current standards and possibilities. *Eur Surg* 40, 19-26.


Katz, A. R., Mukherjee, D. P., Kaganov, A. L. & Gordon, S. (1985). A new synthetic monofilament absorbable suture made from polytrimethylene carbonate. *Surg Gynecol Obstet*, 161, 213-22.

Keck, M., Selig, H. F., Lumenta, D. B., Kamolz, L. P., Mittlbock, M. & Frey, M. (2012). The use of Suprathel((R)) in deep dermal burns: first results of a prospective study. *Burns*, 38, 388-95.

Kumar, S., Ashe, H. A., Parnell, L. N., Fernando, D. J., Tsigos, C., Young, R. J., Ward, J. D. & Boulton, A. J. (1994). The prevalence of foot ulceration and its correlates in type 2 diabetic patients: a population-based study. *Diabet Med*, 11, 480-4.

Larson, M.L., Elkady, D., Sharma, S., Beaucock, B., Lou, R.B., Khandelwal, A. (2024). *Burns*. 2024 Sep;50(7):1832-1839

Liden BA, Ramirez-GarciaLuna JL. Efficacy of a polylactic acid matrix for the closure of Wagner grade 1 and 2 diabetic foot ulcers: a single-center, prospective randomized trial. *Index Wounds*. 2023;35(8):E257-E260

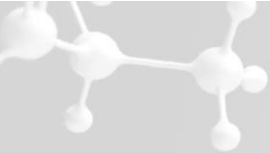Liden BA, Liu L, Regulski M, Foster M, Leon R, Palazzi G, Ramirez-GarciaLuna JR. A multicenter retrospective study comparing a polylactic acid CAMP with intact fish skin graft or a collagen dressing in the management of diabetic foot ulcers and venous leg ulcers. *Wounds*. 2024 Sep;36(9):297-302

Lindford, A. J., Kaartinen, I. S., Virolainen, S. & Vuola, J. (2011). Comparison of Suprathel(R) and allograft skin in the treatment of a severe case of toxic epidermal necrolysis. *Burns*, 37, e67-72.

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März V, Vogt M. Skin Healing of Deep Second Degree Burn Injuries in Four Individuals Sustained in a Boat Explosion Results after Different Approaches. *Eur. Burn J*. 2020, 1, 191–195

Merz, K. M., Sievers, R., Reichert, B. (2011). Suprathel® for coverage of superficial dermal burns of the face. *GMS Verbrennungsmedizin*, 4.

Miguel-Ferrero M, Delgado-Miguel C, Díaz M, Carlos López-Gutiérrez J. Toxic epidermal necrolysis management with suprathel™ Tratamiento de la necrólisis epidérmica tóxica con suprathel®. *An Pediatr (Engl Ed)*. 2023 Aug 17:S2341-2879(23)00185-0


Moellhoff N, Lettner M, Frank K, Giunta RE, Ehrl D. Polylactic Acid Membrane Improves Outcome of Split-Thickness Skin Graft Donor Sites: A Prospective, Comparative, Randomized Study. *Plast Reconstr Surg*. 2022 Nov 1;150(5):1104-1113

Moffatt, C. J. & Dorman, M. C. (1995). Recurrence of leg ulcers within a community ulcer service. *J Wound Care*, 4, 57-61.

Moher D, Liberati A, Tetzlaff J, Altman DG, The PRISMA Group. Preferred Reporting Items for Systematic Reviews and Meta-Analyses: The PRISMA Statement. *PLoS Med*. 2009;6:e1000097. DOI:10.1371/journal.pmed1000097

Mueller, E., Haim, M., Petnehazy, T., Acham-Roschitz, B. & Trop, M. (2010). An innovative local treatment for staphylococcal scalded skin syndrome. *Eur J Clin Microbiol Infect Dis*, 29, 893-7.

Nischwitz SP, Popp D, Shubitidze D, Luze H, Zrim R, Klemm K, Rapp M, Haller HL, Feisst M, Kamolz LP. The successful use of polylactide wound dressings for chronic lower leg wounds: A retrospective analysis. *Int Wound J*. 2021 Nov 8.

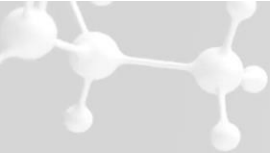O'meara, S. & Martyn-St James, M. (2013). Foam dressings for venous leg ulcers. Cochrane Database Syst Rev, CD009907.

Pego, A. P., Van Luyn, M. J., Brouwer, L. A., Van Wachem, P. B., Poot, A. A., Grijpma, D. W. & Feijen, J. (2003). In vivo behavior of poly(1,3-trimethylene carbonate) and copolymers of 1,3-trimethylene carbonate with D,L-lactide or epsilon-caprolactone: Degradation and tissue response. J Biomed Mater Res A, 67, 1044-54.

Pfurtscheller, K. & Trop, M. (2014). Phototoxic plant burns: report of a case and review of topical wound treatment in children. Pediatr Dermatol, 31, e156-9.

Pfurtscheller, K., Zobel, G., Roedl, S. & Trop, M. (2008). Use of Suprathel dressing in a young infant with TEN. Pediatr Dermatol, 25, 541-3.

Pitt, C. G., Gratzl, M. M., Kimmel, G. L., Surles, J. & Schindler, A. (1981). Aliphatic polyesters II. The degradation of poly (DL-lactide), poly (epsilon-caprolactone), and their copolymers in vivo. Biomaterials, 2, 215-20.

Quinn, K. J., Courtney, J. M., Evans, J. H., Gaylor, J. D. & Reid, W. H. (1985). Principles of burn dressings. Biomaterials, 6, 369-77.

Rahmanian-Schwarz, A., Beiderwieden, A., Willkomm, L.M., Amr, A., Schaller, H.E., Lotter, O. A clinical evaluation of Biobrane® and Suprathel® in acute burns and reconstructive surgery. Burns. 2011 Dec;37(8):1343-8

Rajendran, S., Anand, S.C.. (2011). Hi-tech textiles for interactive wound therapies: Handbook of Medical Textiles.

Rashaan, Z. M., Krijnen, P., Allema, J. H., Vloemans, A. F., Schipper, I. B. & Breederveld, R. S. (2017). Usability and effectiveness of Suprathel® in partial thickness burns in children. Eur J Trauma Emerg Surg, 43, 549-556.

Robson, M. C., Steed, D. L. & Franz, M. G. (2001). Wound healing: biologic features and approaches to maximize healing trajectories. Curr Probl Surg, 38, 72-140.

Rothenberger, J., Constantinescu, M. A., Held, M., Aebbersold, D. M., Stolz, A., Tschumi, C. & Olariu, R. (2016). Use of a Polylactide-based Copolymer as a Temporary Skin Substitute for a Patient With Moist Desquamation Due to Radiation. Wounds, 28, E26-30.

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
Ruckley, C. V. (1998). Caring for patients with chronic leg ulcer. BMJ, 316, 407-8.

Sari, E., Eryilmaz, T., Tetik, G., Ozakpinar, H. R. & Eker, E. (2014). Suprathel® -assisted surgical treatment of the hand in a dystrophic epidermolysis bullosa patient. Int Wound J, 11, 472-5.

Schiefer, J.L., Rahmanian-Schwarz, A., Schaller, H.E., Manoli, T. A Novel Hand-shaped Suprathel simplifies the Treatment of Partial-Thickness Burns. Adv Skin Wound Care. 2014 Nov;27(11):513-6

Schiefer JL, Andreae J, Bagheri M, Fuchs PC, Lefering R, Heitzmann W, Schulz A. A clinical comparison of pure knitted silk and a complex synthetic skin substitute for the treatment of partial thickness burns. Int Wound J. 2022 Jan;19(1):178-187.

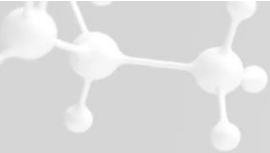Schiefer JL, Aretz GF, Fuchs PC, Bagheri M, Funk M, Schulz A, Daniels M. Comparison of wound healing and patient comfort in partial-thickness burn wounds treated with SUPRATHEL and epicte hydro wound dressings. *Int Wound J*. 2022 May;19(4):782-790.

Schiefer JL, Aretz FG, Fuchs PC, Lefering R, Yary P, Opländer C, Schulz A, Daniels M. Comparison of Long-Term Skin Quality and Scar Formation in Partial-Thickness Burn Wounds Treated with Suprathel® and epicitehydro® Wound Dressings. *Medicina (Kaunas)*. 2022 Oct 28;58(11):1550. (b)

Schiefer JL, Andreae J, Fuchs PC, Lefering R, Heidekrueger PI, Schulz A, Bagheri M. Evaluation of Scar Quality after Treatment of Superficial Burns with Dressilk® and Suprathel®-In an Intraindividual Clinical Setting. *J Clin Med*. 2022 May 18;11(10):2857 (c)

Schreml, S., Szeimies, R. M., Prantl, L., Karrer, S., Landthaler, M. & Babilas, P. (2010). Oxygen in acute and chronic wound healing. *Br J Dermatol*, 163, 257-68.

Schriek K, Ott H, Sinnig M. Paradigm Shift in Treatment Strategies for Second-Degree. *Eur. Burn J*. 2022, 3, 1–9

Schwarze, H., Kuntscher, M., Uhlig, C., Hierlemann, H., Prantl, L., Noack, N. & Hartmann, B. (2007). Suprathel, a new skin substitute, in the management of donor sites of split-thickness skin grafts: results of a clinical study. *Burns*, 33, 850-4.

Schwarze, H., Kuntscher, M., Uhlig, C., Hierlemann, H., Prantl, L., Ottomann, C. & Hartmann, B. (2008). Suprathel, a new skin substitute, in the management of partial-thickness burn wounds: results of a clinical study. *Ann Plast Surg*, 60, 181-5.

Selig, H. F., Keck, M., Lumenta, D. B., Mittlbock, M. & Kamolz, L. P. (2013). The use of a polylactide-based copolymer as a temporary skin substitute in deep dermal burns: 1-year follow-up results of a prospective clinical noninferiority trial. *Wound Repair Regen*, 21, 402-9.

Sun, H., Mei, L., Song, C., Cui, X. & Wang, P. (2006). The in vivo degradation, absorption and excretion of PCL-based implant. *Biomaterials*, 27, 1735-40.

Szycher, M. & Lee, S. J. (1992). Modern wound dressings: a systematic approach to wound healing. *J Biomater Appl*, 7, 142-213.

Tams, J., Joziassse, C. A., Bos, R. R., Rozema, F. R., Grijpma, D. W. & Pennings, A. J. (1995). High-impact poly(L/D-lactide) for fracture fixation: in vitro degradation and animal pilot study. *Biomaterials*, 16, 1409-15.


Thomas, S. S., Lawrence, J. C. & Thomas, A. (1995). Evaluation of hydrocolloids and topical medication in minor burns. *J Wound Care*, 4, 218-20.

Trabold, O., Wagner, S., Wicke, C., Scheuenstuhl, H., Hussain, Z., Rosen, N., Seremetiev, A., Becker, H., Hunt, T (2003). Lactate and oxygen constitute a fundamental regulatory mechanism in wound healing. *Wound Rep Reg* 11:504–509.

Uhlig, C., Hierlemann, H., Dittel, K.-K. (2007). Actual Strategies in the Treatment of Severe Burns - Considering Modern Skin Substitutes. *Osteo trauma care* 15, 2-7.

Uhlig, C., Rapp, M. & Dittel, K. K. (2007a). [New strategies for the treatment of thermally injured hands with regard to the epithelial substitute Suprathel]. *Handchir Mikrochir Plast Chir*, 39, 314-9.

Uhlig, C., Rapp, M., Hartmann, B., Hierlemann, H., Planck, H. & Dittel, K. K. (2007b). Suprathel-an innovative, resorbable skin substitute for the treatment of burn victims. *Burns*, 33, 221-9.

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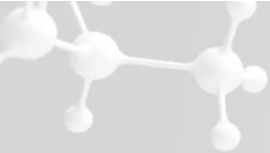Vowden, K. R. & Vowden, P. (2009). The prevalence, management, equipment provision and outcome for patients with pressure ulceration identified in a wound care survey within one English health care district. *J Tissue Viability*, 18, 20-6.


Wasiak, J., Cleland, H., Campbell, F. & Spinks, A. (2013). Dressings for superficial and partial thickness burns. *Cochrane Database Syst Rev*, CD002106.

Wasiak, J., Cleland, H. (2005). Minor thermal burns. *Clinical Evidence*, 2754-2763.

Wild, S., Roglic, G., Green, A., Sicree, R. & King, H. (2004). Global prevalence of diabetes: estimates for the year 2000 and projections for 2030. *Diabetes Care*, 27, 1047-53.

Woodward, S. C., Brewer, P. S., Moatamed, F., Schindler, A. & Pitt, C. G. (1985). The intracellular degradation of poly(epsilon-caprolactone). *J Biomed Mater Res*, 19, 437-44.

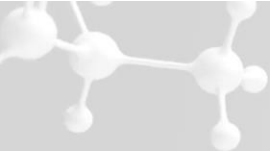Wu, L., Norman, G., Dumville, J. C., O'meara, S. & Bell-Syer, S. E. (2015). Dressings for treating foot ulcers in people with diabetes: an overview of systematic reviews. *Cochrane Database Syst Rev*, CD010471.

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